





# Meditech Europe

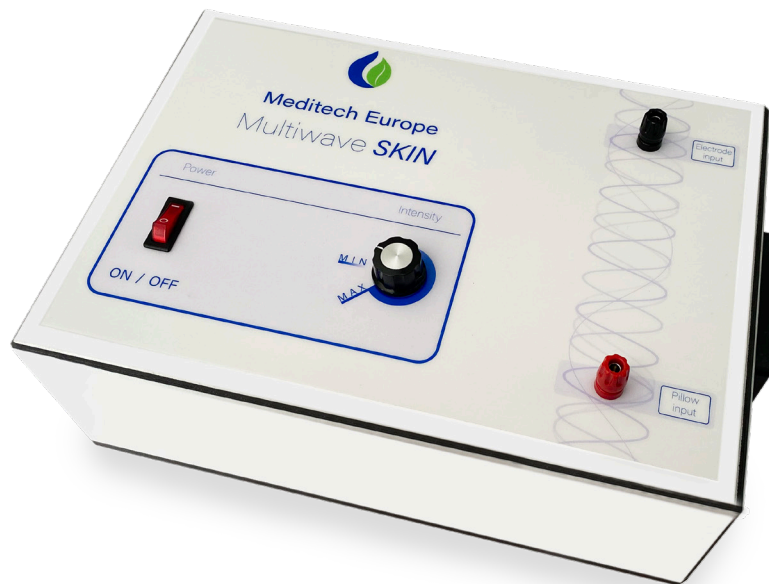

Multiwave Skin usermanual

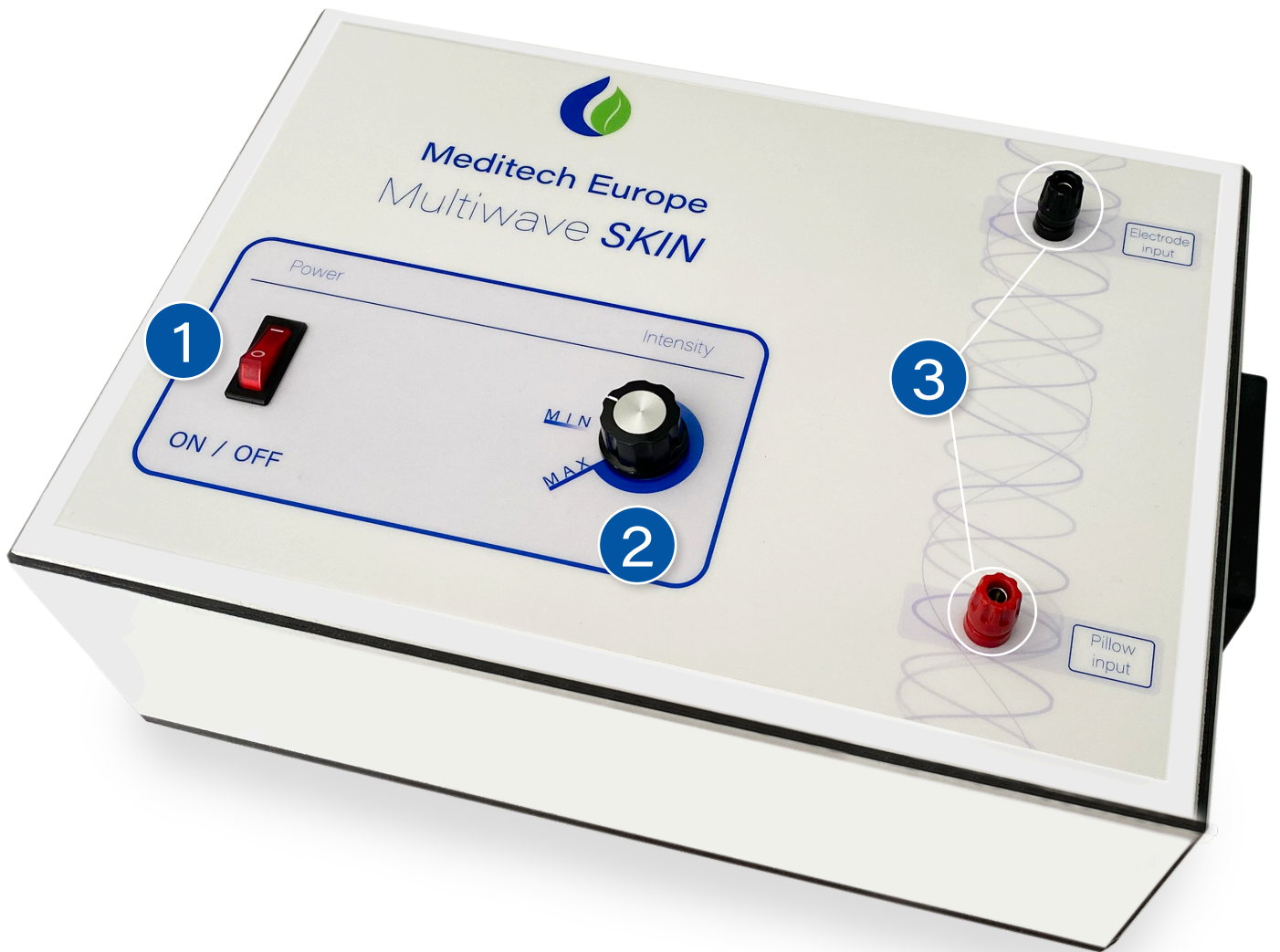

# Operating the Multiwave SKIN

1

## **On / Off button**

Switches the Multiwave Skin on and off.

On: The light in the button is on.

Off: The light in the button will not illuminate.

2

## **Intensity**

Adjusts the Multiwave Skin strength.

When the button is set to Minimal (Min), the strength is low. When the knob is at Maximum (Max) the strength is maximum.

We advise you to start every treatment with the setting at Min and slowly increase it until the used accessory reacts properly. You can read more about this in the explanation of the accessories.

3

## **Connectors**

To connect the accessories to the Multiwave Skin.

The connection for the therapy pillow (Pillow input) is there for the grounding of the user. When the therapy pillow is connected the user takes place on the pillow. Hereby the user becomes a conductor of the energy created by the Multiwave Skin.

The connection for the accessories (Electrode input) is for the connection of the handle with which the treatment of the user is carried out. The accessories are connected to the handle (electrodes) or the Rife Therapy Tube is connected directly to the Multiwave Skin.





## Accessories




The **Rife Therapy Tubes** are very powerful treatment lamps that with their extreme transfer of frequencies and warmth provide a very soothing feeling. When these tubes are connected to the Multiwave Oscillator it is important to also use the Therapy Cushion. It is also important that the setting of the Multiwave Skin is slowly increased until the Rife Therapy Tube is experienced as comfortable. Use the Rife Therapy Tube as often as you wish.

- Connect the Rife Therapy Tube to the Multiwave Skin by plugging the cord into the red connector. Connect the Therapy Cushion to the black connector.
- Sit on the cushion and turn on the Multiwave Skin. Initially set it as low as possible. The Rife Therapy Tube will now turn a light orange.
- Now place the Rife Therapy Tube on the skin in the area to be treated. Where the Rife Therapy Tube touches the skin the colour will be more intense orange and a gentle crackling will be heard. There may also be some ozone smell. This is not harmful!
- Now the Rife Therapy Tube works and transmits the frequencies to your skin, nerve system and blood stream.
- You now move slowly and relaxed counter-clockwise over the skin. If it feels more comfortable you can set the Multiwave a little higher, but make sure that the feeling remains pleasant, without stinging.
- You treat for 15 minutes. If there are several areas that require attention, you can also alternate, but not too often and not too quickly. Remember that the warmth and massaging effect is part of the treatment.

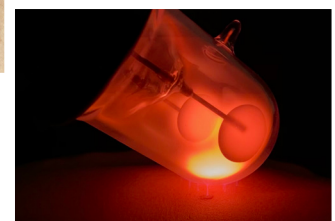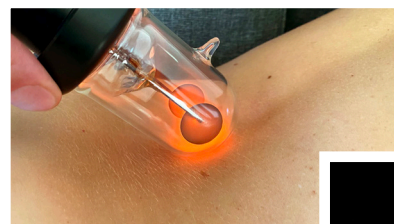



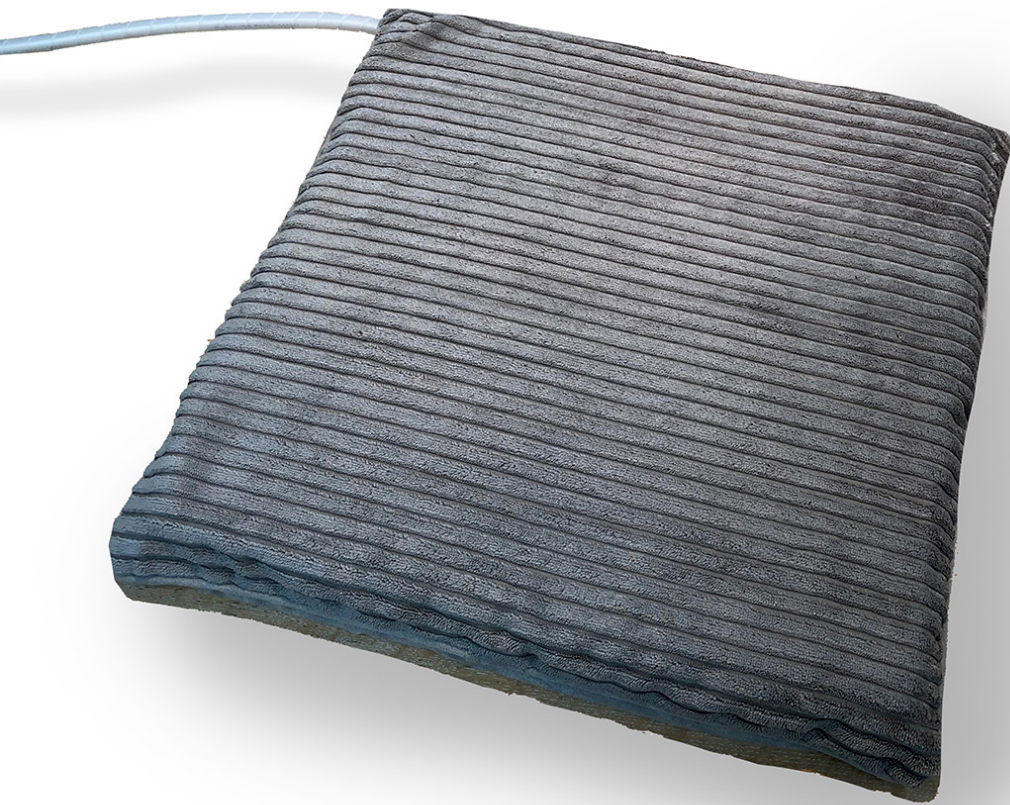

With the **Therapy cushions** you can perform an additional treatment. You connect both cushions to the low setting Multiwave Skin and take place on one cushion. With the other cushion you can, for example, treat your belly by placing the second cushion over your belly or your feet by placing your feet on the second cushion. The additional Therapy Cushion treatment can be done for 2 x 30 minutes.

In addition you use the Therapy Cushion to ground the treatment accessories, such as the electrodes and the Rife Therapy Tube. Connect the therapy pillow to the Pillow input connector. Connect the other connector to the accessory of your choice. Then take your place on the pillow and switch on the Multiwave Skin. With the Intensity button you can set the Multiwave Skin to the correct strength so that the treatment can proceed optimally.



# Electrodes for the Multiwave Skin

## Multiwave Skin Glasselectrodes

Meditech Europe has a wide range of glass electrodes that can be used in combination with the Multiwave Skin.

To work with the electrodes it is important to know how to use them.

## Connecting

To work with the glass electrodes, you need a hand unit, to which the glass electrodes are connected. When connecting the glass electrodes, make sure that this is done accurately and carefully due to the fragility of the glass electrodes.

When connecting the cable from the hand unit to the Multiwave Skin (Electrode input) make sure that the cushion is connected to the input (Pillow input) and that the user actually sits on the cushion. This in connection with the grounding. If you do not do this, the electrode cannot work.



When you are ready to use the electrode, switch on the Multiwave Skin.

You have the hand unit with the glass electrode in your hand and turn very calmly and slowly the Intensity knob slightly open, until you hear a very soft sparking sound.

When you hold your hand against the glass electrode you will see a light glow in the electrode and you will feel a slight tingling. This should be a pleasant tingling and in no way an unpleasant experience. If there is an unpleasant feeling, you should adjust the Intensity knob until a comfortable position is reached.

With this setting, you can now begin treatment of the affected area. You can continue this treatment until the user indicates that it is sufficient. Maximum 10-15 minutes treatment on one spot.

To make the treatment more pleasant, you can use a lotion to make the conduction more pleasant. Make sure that this lotion does not contain any alcohol. A good product is Colloidal Ointment or S&G BodyLotion from Meditech Europe. This is a 100% natural lotion that you can use for this purpose.



**Electrode 0**

Flat electrode with larger surface area for abdomen, stomach



**Electrode 1a**

Small flat electrode for gentle facial massage



**Electrode 2t**

Amplifier electrode, special for the treatment of rheumatism



**Electrode 2t groot**

Amplifier electrode, with legs, for extra powerful operation



**Electrode 3**

Neck and neck electrode for external treatment



**Electrode 5**

Comb electrode for hair care and scalp treatment



**Electrode 6**

Spine electrode, isolated



**Electrode 7**

Ooglideelectrode



**Electrode 8**

Eyelid electrode



**Electrode 9**

Eye electrode, double shape



# Electrodes voor de Multiwave Skin



**Electrode 10**  
Fulguration electrode  
for the removal of warts,  
isolated



**Electrode 11**  
Nose electrode, isolated



**Electrode 12**  
Ear electrode, Insulated



**Electrode 13**  
Throat electrode, isolated



**Electrode 14**  
Glass-saturator or satiator



**Electrode 15**  
Tongue electrode, isolated



**Electrode 16**  
Nerv electrode



**Electrode 17**  
Urethra electrode,  
insulated



**Electrode 18**  
Penis electrode, fork  
shape



**Electrode 20**  
Vagina electrode, isolated



# Electrodes voor de Multiwave Skin



**Electrode 24**  
Endothermal electrode,  
isolated



**Electrode 25**  
Prostate electrode,  
isolated



**Electrode 26**  
Isolated prostate electrode



**Electrode 30**  
Spiral shaped heart  
electrode



**Electrode 31a**  
Electrode for the  
development of the breast,  
with bulb



**Electrode 31b**  
Electrode for the body,  
with bulb



**Electrode 37**  
Tooth electrode with  
spherical cup, especially  
for the treatment of  
pyorrhoea and Gingivitis



**Electrode 39**  
Pointed tooth electrode for  
treatment of abscesses



**Electrode 51n**  
Cold neon massage roller,  
larger areas



**Electrode 101**  
Ozone inhaler with cotton  
ball





# Meditech Europe BV

Daalder 14

8305 BE Emmeloord

tel.: +31 527 292331

mail: [info@meditecheurope.nl](mailto:info@meditecheurope.nl)

web: [www.meditecheurope.nl](http://www.meditecheurope.nl)