Meditech Europe

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fig. 2

Package Contents

- Soladey Ion 5 handle
- Replaceable brush head (1 x Honeycomb)
- Instruction manual

Before Starting

- Carefully read the instruction manual.
- Ensure all parts are present.
- Wash the brush head and handle with water before first use.

Getting Started

1. Attaching the Brush Head: Firmly attach a brush head to the Soladey Ion 5 handle. Ensure it is securely in place.
2. Preparation: Wet the brush head with water. This activates the ionic cleaning technology of the Soladey Ion 5.
3. Optional: You may use toothpaste as desired, although the ionic technology can also clean effectively without toothpaste.

Brushing Your Teeth

- Brush: Move the toothbrush in gentle, circular motions along all surfaces of your teeth and molars. Pay special attention to the back of your teeth and hard-to-reach places.
- Duration: Brush for at least two minutes. Some models feature a built-in timer to assist with this.
- After Brushing: Rinse your mouth and the brush head thoroughly with water.

Maintenance and Care

- After Each Use: Rinse the brush head thoroughly under running water.
- Storage: Store the toothbrush in a clean, dry place, preferably standing upright so the brush head can dry.
- Replacing the Brush Head: Replace the brush head every 3 months, or sooner if the bristles begin to spread apart or show signs of wear.
- For all variants of replaceable brush heads, check out Meditech Europe's website.

Safety Information

- Soladey5 is not an electrical toothbrush.
- People who use an Electric medical device, etc. should not use this product.
- Keep the brush head attached to the handle and rod at all times except during replacement.
- Never leave the toothbrush or bare handle and rod near young children or infants.
- Keep the brush head attached to the handle and rod at all times except during replacement.
- Never leave the toothbrush or bare handle and rod near young children or infants.

Troubleshooting

- **Toothbrush Not Working:** Ensure the brush head is correctly attached and has been wetted.
- **Insufficient Cleaning:** Check if the brush head needs replacing or try a different brushing technique.

Contact and Support

For questions or concerns, please contact Meditech Europe through the website or customer service. Keep your proof of purchase and product details for any warranty claims.

fig. 1

fig. 2

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fig. 1

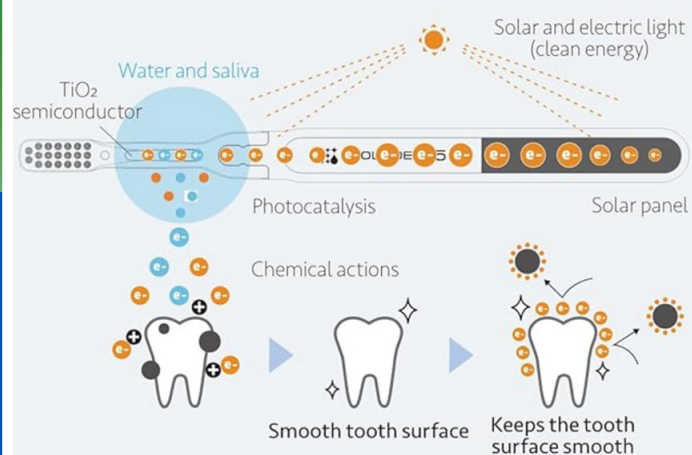
fig. 2

Mechanisms for removing plaque and fighting off cavity-causing bacteria

 Meditech Europe

Meditech Europe

Daalder 14, 8305 BE Emmeloord
tel.: +31 527 292331
mail: info@meditecheurope.nl
web: www.meditecheurope.nl