Meditech Europe

Soladey user manual

Meditech Europe

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Meditech Europe

The SOLADEY ION3 uses the same technology and principles as the SOLADEY-Eco with the added effectiveness of a solar panel at the end of the handle. This makes it more effective in catching light and therefore able to function in low-light situations and with the mouth closed.

Instructions for use

1. Before using the SOLADEY toothbrush, moisten the brush and semiconductor with water. For outdoor use saliva is also sufficient*.
2. First brush the back (not the throat!) and then the front of the tongue, without using toothpaste. Especially in the morning, this will remove the stomach acids you have accumulated while sleeping. You may notice discolouration on the bristles (concentrated bacteria and accumulation of acids). To do this, rinse the brush repeatedly with water. The concentration of bacteria in your mouth is now lower. Please note that the tongue generally contains the majority of oral bacteria! If you consider the brushes and tongue sufficiently clean, you can work on the tooth surfaces.
3. Brush your teeth lightly, in the same way as you would with a regular toothbrush. Concentrate on the contact with the teeth without using too much force. For this reason we recommend a MEDIUM toothbrush. In contrast to a regular toothbrush, you do not need to brush by force. The medium brush also brushes in the spaces between your teeth and works more efficiently than brushing with a soft head. Brushing lightly is better for your gums and helps fight gingivitis (gum disease / receding gums).

*Note: Saliva is sufficient for brushing teeth outside the home. We do not recommend brushing outside the home if there is no source of running water available for cleaning the brushes.

An adequate light source (natural or artificial) is required to activate the semiconductor. If you want to brush the surfaces at the back of the teeth, there is limited direct light. If you keep your mouth closed, thus covering the ionic rod, there is a lack of light and no ionic reaction takes place. Be aware that when you brush the back of the teeth, the bar must receive enough light. Open your mouth as wide as possible, so that the available light can reach the back of your teeth. The best option is to brush your teeth in the direct vicinity of a light source. That light source will illuminate the inside of your mouth. The stronger the light source, the better the result. When you use the toothbrush effectively, you will notice that your saliva will foam slightly. If you use toothpaste, use very little. This can restrict the brush from reaching lightly.

Be sure to clean the toothbrush after each use. If you live in an area with hard water, carefully remove the brush from the handle. You can clean the handle with brushes (use your old toothbrush or another head).


Meditech Europe

We recommend using natural mouthwash and toothpaste as a substitute for freshening up your breath. A few drops of NATURAL MINT OIL or even OREGANO OIL also works well.

SAFETY: Make sure the brush is attached to the handle and to the rod (except during replacement). Keep the toothbrush, the bare handle and/or rod away from children!
Replacement of the brush

In order for your SOLADEY toothbrush to work as efficiently as possible, you should replace the brush regularly. We recommend changing the brush every 2 to 3 months or when the bristles are standing out. The figure below shows the correct method of changing the brush:

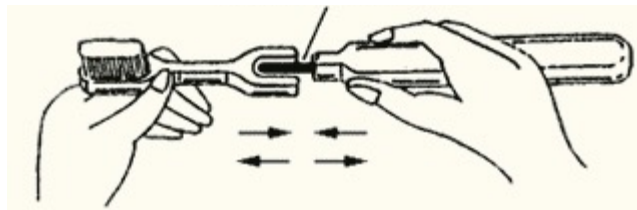Detach: Pull the brush (head) from the handle, along the titanium rod.
Attach: Insert the head into the handle.

DO NOT twist and/or bend the head of the toothbrush when detaching and attaching the brush.

Recommendations

- Rinse the toothpaste completely off the toothbrush after use.
- Do not remove the brush from the handle except when replacing it.
- Do not leave it in the mouth without a brush head.